

## Directions for A-D (Below Knee) Garments

Videos may be available on our website. You must create an account and login to view them. **Live web-based in-services are also available. We highly recommend that you schedule an appointment for an in-service.**

A practitioner's order is required for all garments. DME stores do not need to forward a practitioner's order to Farrow Medical. Have the patient sign the authorization. Please note that OTS sizes can accommodate at least a 15 to 20% reduction in edema, but only a 5% increase.

Fill in the Patient / Billing & Shipping Information areas **completely**. Next, choose the FarrowWrap™ version - Classic or LITE. Classic is generally reserved for more severe edema cases and is only available in tan. LITE is more appropriate for mild to moderate edema and is available in tan and black. On the measurement page, part of the sku number (ie FW\_-O-LR) is represented by a "\_" which represents either a "CL" for Classic or "LT" for LITE. The appropriate sku number can be cross-referenced with the Pricing Schedule to calculate the price.

Only circumferences at points **A**, **B**, & **C** and length measurements **X** & **A-D** are necessary to determine if the patient can fit into an **OTS (Off-The-Shelf or ready-made) A-D garment**. If the patient will not fit into an OTS garment, then the additional measurements for that particular component are required. All length measurements are measured along the posterior contour of the leg following contours around any lobes. If any lobes are present, drawings and/or pictures (with the patient's permission) will help. Measure only portions of the limb to be fitted with a garment. Follow these instructions on how to properly measure. **ONLY** record **ACTUAL** measurements.

### FarrowWrap™ Legpiece

- Obtain the following 3 measurements:
  - Mid or widest calf circumference (at **C**).
  - Least ankle circumference (at **B**), just above malleoli or bumps on either side of ankle.
  - Length of leg following posterior contour from 2 finger widths below posterior knee crease to the ground, **A-D**.
- Now, match the OTS midcalf and ankle circumferences to a respective size. If the mid-calf **AND** ankle circumferences do not both match a single size, or if the patient has a large lobule, then a custom legpiece is highly recommended and we will need all associated measurements.
- If the circumferences do match an OTS size, then compare **A-D** to the OTS regular length and tall length measurements. If **A-D** falls within 1 - 1.5 centimeters of one of these lengths, then in the space next to the appropriate OTS length enter the number of legpieces you would like in that given size. If the measured length is not within 1.5 centimeters of either length, a custom legpiece is highly recommended and all associated measurements would be required. Enter the number of legpieces desired.

### Note

The leg measurements used by this form allow for at least a 20% additional reduction in leg and ankle circumference (edema reduction), but can only accommodate about a 5% increase. There is room for greater variance in ankle circumferences than for calf circumferences. In general, Classic will accommodate greater reductions in edema and LITE will accommodate greater increases. Patients may fit into an OTS garment even though this form may indicate otherwise, but we cannot guarantee it.

### FarrowWrap™ Footpiece or FarrowHybrid AD Liner

- Measure the affected foot's midfoot circumference, **A**, and length, **X**. The length is measured from the proximal (posterior) border of either the 1st or 5th metatarsal head (bunion), which ever is shorter, to the posterior most aspect of the heel. This measurement should be made across the floor - not along the side of the patient's foot.
- Next, match **A** to a respective-sized OTS footpiece. If a match cannot be made, a custom footpiece will be required.
- Next, compare **X** to the OTS regular length and long length measurements. If **X** falls within a 0.5 - 1.0 cm of one of these OTS lengths, then in the space next to the appropriate length enter the number of footpieces you would like in that given size. If the length measurement is not within 1.0 centimeter of any OTS length, a custom footpiece is highly recommended.
  - Note: The **X** measurement defaults to extend to the metatarsal heads, for safety reasons. If you wish the footpiece to extend to the base of the toes, then you must add this extra distance to **X** when choosing an OTS garment. For custom footpieces, under "Extend footpiece length to base of toes," write L (left), R (right), or BL (bilateral) on the form and we will add 2 cm to length.  
**Only do this on patients with adequate circulation and no peripheral neuropathy; do not use on symptomatic diabetics.**
- Enter the number of footpieces desired.
- As a replacement option in mild to moderate edema, a FarrowHybrid™ AD Liner may be used in place of a footpiece.

### Liner Options

If no options are selected, a pair of **Farrow Silver AD Liners** will be sent **standard** (free) with each footpiece / legpiece combination **only**. The FLSL-AD2 can accommodate leg circumferences up to 58 cm and the FLSL-AD4 can accommodate leg circumferences from 53-70 cm; an appropriate pair will be issued standard. If the leg circumference is greater than 70 cm, a pair of **TG® Soft AD Liners** will be issued **standard** instead, also free. TG® Soft Liners come in 3 circumferences, S < 40 cm, M 40-70 cm, & L 70-125 cm. They may be **substituted** for silver liners in the case of silver incompatibilities without charge. To order a **different than standard** liner, **check the "Use Nonstandard Liner," and use the liner order form**. To order additional liners, also use the liner order form. The **FarrowFoam™ AD Liner** is for patients with fragile skin for use under the Classic leg and footpieces. They come in 4 sizes and 2 models. See catalog and liner order form for details. Please note **when ordering a FarrowFoam™ Liner for use under an OTS garment, that you must also add 10 cm to leg circumferences and 5 cm to foot circumferences when determining the correct sized FarrowWrap™ leg or foot piece**. We recommend at least 2 pairs of liners per extremity.

### Other Options


Some patients have very triangular shaped legs and even the Classic garment may have difficulty staying up. Additional velcro-like spines, called Hook Stays, can be easily attached to the garment to help hold it in place. The stays are available in small (16 cm length), medium (21 cm length), and large (27 cm length). For custom garments, additional velcro tabs may be ordered to help prevent gapping between bands on oddly shaped limbs.



# FarrowWrap™ AG OTS and Custom Order Form



Business: \_\_\_\_\_ Patient: \_\_\_\_\_ PO: \_\_\_\_\_  
 Person Measuring: \_\_\_\_\_ Phone/Fax: \_\_\_\_\_ / \_\_\_\_\_ Date: \_\_\_\_\_

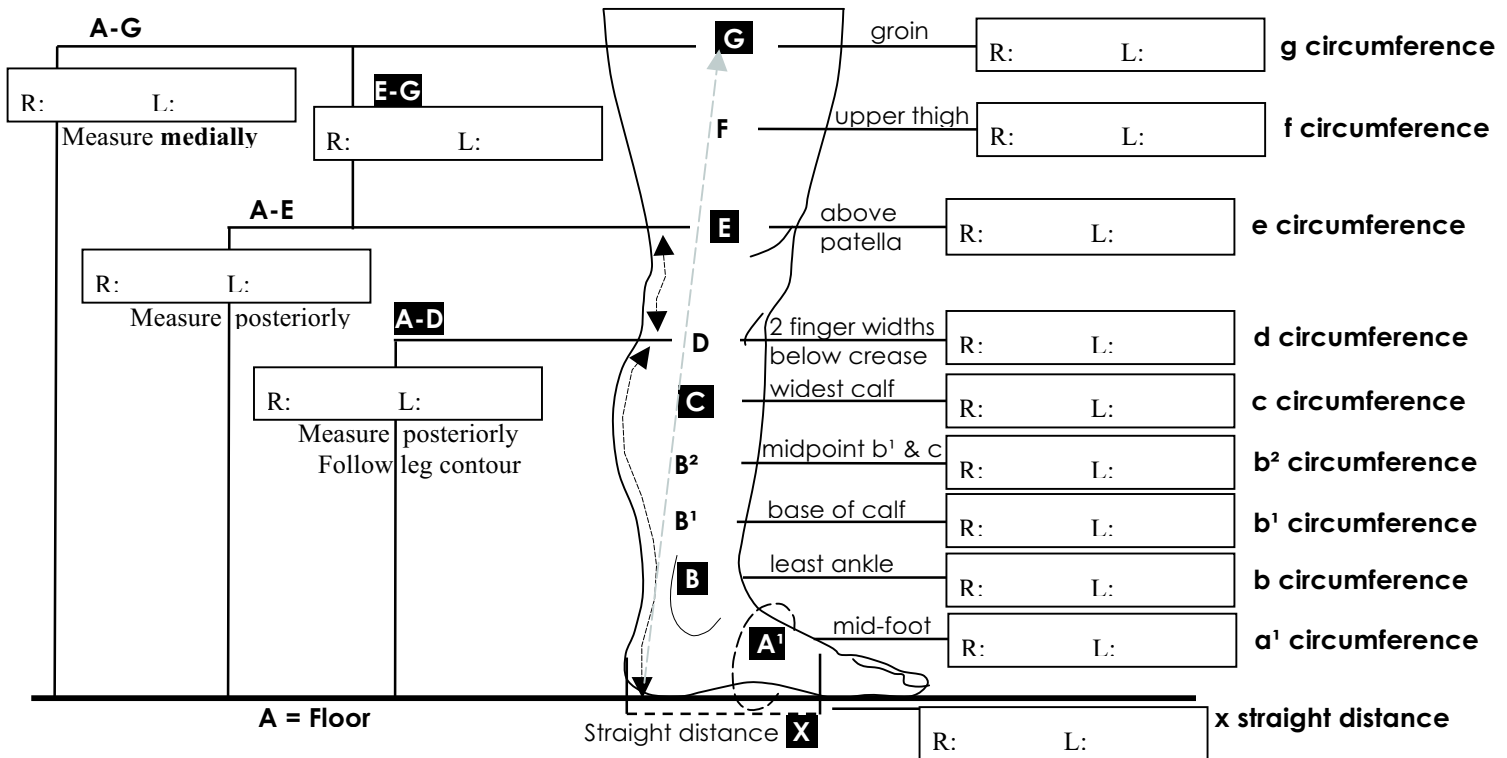

OTS Components		Thighpiece (Only comes with wrap-around kneepiece)									
		XSmall	Qty	Small	Qty	Medium	Qty	Large	Qty	XLarge	Qty
(G) Groin		60–70 cm		65–75 cm		70–80 cm		75–90 cm		85–100 cm	
(E) Above Knee		47–52 cm		53–59 cm		60–66 cm		67–73 cm		74–80 cm	
(E-G) Short	FW_-O-TS	18 cm		18 cm		18 cm		18 cm		18 cm	
(E-G) Regular	FW_-O-TR	23 cm		23 cm		23 cm		23 cm		23 cm	
(E-G) Tall	FW_-O-TT	29 cm		29 cm		29 cm		29 cm		29 cm	
		Legpiece									
		XSmall	Qty	Small	Qty	Medium	Qty	Large	Qty	XLarge	Qty
(C) Mid-Calf		36–43 cm		42–50 cm		48–58 cm		53–63 cm		58–68 cm	
(B) Ankle		21–25 cm		25–30 cm		30–36 cm		36–42 cm		42–50 cm	
(A-D) Regular	FW_-O-LR	35 cm		37 cm		39 cm		41 cm		41 cm	
(B-D)		28 cm		30 cm		31 cm		33 cm		33 cm	
(A-D) Tall	FW_-O-LT	40 cm		42 cm		44 cm		46 cm		46 cm	
(B-D)		33 cm		35 cm		36 cm		38 cm		38 cm	
		Footpiece									
		XSmall	Qty	Small	Qty	Medium	Qty	Large	Qty	XLarge	Qty
(A') Mid-Foot		22–24 cm		25–27 cm		28–30 cm		31–34 cm		35–40 cm	
(X) Regular	FW_-O-FR	16.5 cm		18 cm		19.5 cm		21 cm		23 cm	
(X) Long	FW_-O-FL	18.5 cm		20 cm		21.5 cm		23 cm		25 cm	

Options	
Item	√
Classic (Tan)	
LITE (Tan)	
Use Nonstandard Liner Option (see Liner Form)	
When a footpiece & legpiece are purchased together, 1 pair of <b>standard</b> liners are included automatically free. Other liners may cost extra. See directions, liner form, & pricing for details.	

Custom Components	Thighpiece		<input type="checkbox"/> Wrap-around <input type="checkbox"/> Standard		Legpiece		Footpiece		Lymphedema Footpiece Option	
	Right:	Left:	Right:	Left:	Right:	Left:	Right:	Left:	L R BL:	
Quantity										

Other: \_\_\_\_\_